

0.6–6.0 mg/L NO_2^- -N or
2.0–20.0 mg/L NO_2^- High Range

TNTplus®—Method 10237

Scope and application: For wastewater, drinking water, surface water and process water.



Test preparation

Reagent storage

Storage temperature: 15–25 °C (59–77 °F)

pH/Temperature

The pH of the water sample must be between pH 3–10.

The temperature of the water sample and reagents must be between 15–25 °C (59–77 °F).

Before starting

Not more than **3 hours** should elapse between sampling and analyzing the sample.

Review safety information and expiration date on the package.


Review the Safety Data Sheets (MSDS/SDS) for the chemicals that are used. Use the recommended personal protective equipment.

Dispose of reacted solutions according to local, state and federal regulations. Refer to the Safety Data Sheets for disposal information for unused reagents. Refer to the environmental, health and safety staff for your facility and/or local regulatory agencies for further disposal information.


Procedure




1. Carefully remove the foil from the screwed-on DosiCap Zip.

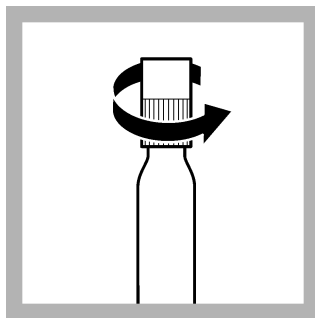


2. Unscrew the DosiCap Zip.

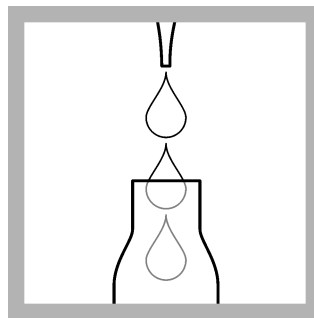


3. Carefully pipet **0.2 mL** of **sample**.

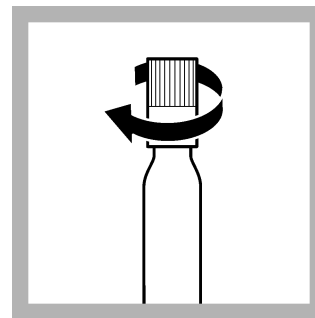


4. Immediately screw the DosiCap Zip back on; **fluting at the top**.

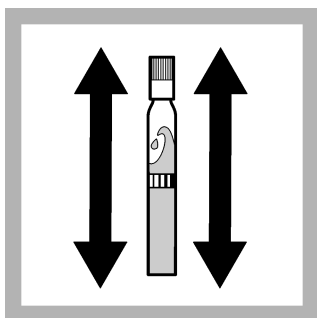

5. Shake **vigorously** until the freeze-dried contents are **completely dissolved**.

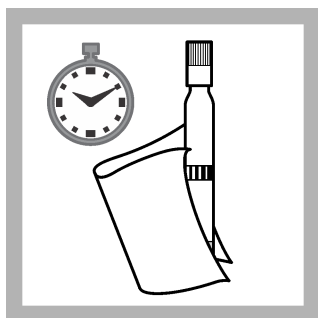

6. After **10 minutes**, thoroughly clean the outside of the vial and evaluate.

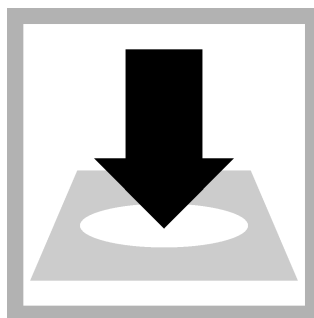

7. Insert the vial into the cell holder.
DR 1900: Go to LCK/TNTplus methods. Select the test, push **READ**.

Interferences

The ions listed in the table have been individually checked against the given concentrations and do not cause interference. The cumulative effects and the influence of other ions have not been determined. Chromium(VI) ions interfere with the determination. Copper(II) ions interfere with the determination even at concentrations below 1 mg/L.

The measurement results must be subjected to plausibility checks (dilute and/or spike the sample).

Interference level	Interfering substance
4000 mg/L	SO_4^{2-}
2000 mg/L	K^+ , NO_3^- , Ca^{2+} , Cl^-
1000 mg/L	NH_4^+ , PO_4^{3-}
200 mg/L	Mg^{2+}
100 mg/L	Cr^{3+} , Hg^{2+}
50 mg/L	Co^{2+} , Zn^{2+} , Cd^{2+} , Mn^{2+}
20 mg/L	Fe^{3+} , Ni^{2+} , Ag^+ , Fe^{2+}
10 mg/L	Sn^{4+}

Summary of method

Nitrites react with primary aromatic amines in acidic solution to form diazonium salts. These combine with aromatic compounds that contain an amino group or a hydroxyl group to form intensively colored azo dyes.

TNT^{plus}



FOR TECHNICAL ASSISTANCE, PRICE INFORMATION AND ORDERING:
In the U.S.A. – Call toll-free 800-227-4224
Outside the U.S.A. – Contact the HACH office or distributor serving you.
On the Worldwide Web – www.hach.com; E-mail – techhelp@hach.com

HACH COMPANY
WORLD HEADQUARTERS
Telephone: (970) 669-3050
FAX: (970) 669-2932